



How to choose the correct receiver length for Phonak

Grab a ruler and get ready to measure your best receiver length. Make sure you measure from **0** where your ear joins your head.

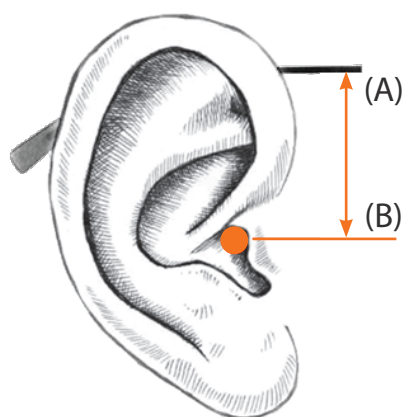

Use a pair of glasses **or a sheet of paper**

Measure the distance from the underside of your glasses frame (A) to the top third of your ear canal (B)

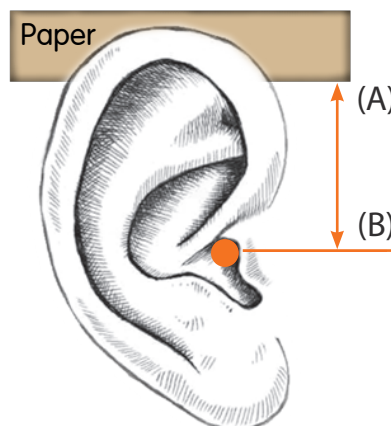

Rest a piece of paper horizontally on top of your ear, and measure from the bottom of the paper (A) to the top third of your ear canal (B)

Phonak hearing aid receiver length guide



SIZE 0	SIZE 1	SIZE 2
 <p>19-23mm</p>	 <p>23-27mm</p>	 <p>27-34mm</p>

If you fit within two ranges, pick the larger size.